



Smart Technology,
Made to Care



**Innovative, Interoperable
Telecare Solutions**

Connect Quick Reference Guide

sensorium.co.uk

Contents

- Page 3 - Switching On/Off
- Page 4 - Programming the Modes/Options
- Page 5 - Modes/Options
- Page 6 - Absence Delay
- Page 7 - Registering to a Notifier / Other System




Smart Technology,
Made to Care

[Connect

Switching On / Off

To switch on, press and hold the external switch. The LED Indicator will flash red & green and then produce a high-pitched tone followed by a brief green colour flash which confirms it is on. Now release the switch.

To switch off, press and hold the external switch. The LED Indicator will flash red & green and then produce a low-pitched tone followed by a brief red colour flash which confirms it is off. Now release the switch.

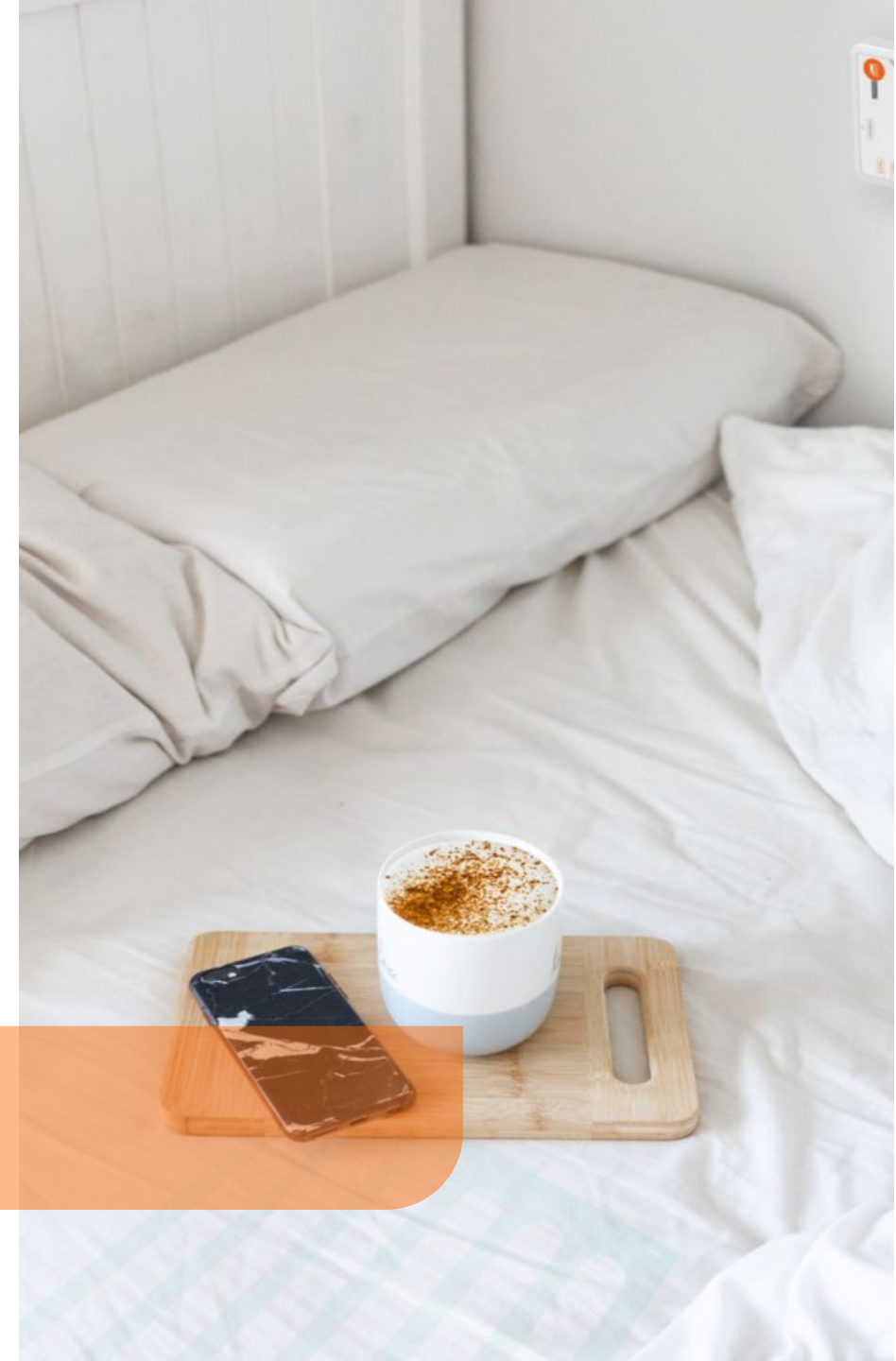


Smart Technology,
Made to Care

Connect

Programming the Modes / Options

To program the modes/options (which can be found on the next page), use the following method:

- Turn the dial to the desired setting
- Press and hold the internal switch using a suitable instrument until LED 2 changes colour and a sound is heard (approx 2 seconds)
- Release the internal switch. LED 2 and the sounder will flash/beep to confirm the mode number you have chosen. Please note there is an extra step if using modes 3 and 4, see the Connect programming guide for more details.
- Repeat steps 1-3 for other modes/options
- If using Modes 1-4, simply set the dial to the desired absence delay as described on the Absence Delay page. There is no need to press the internal switch to set the absence delay

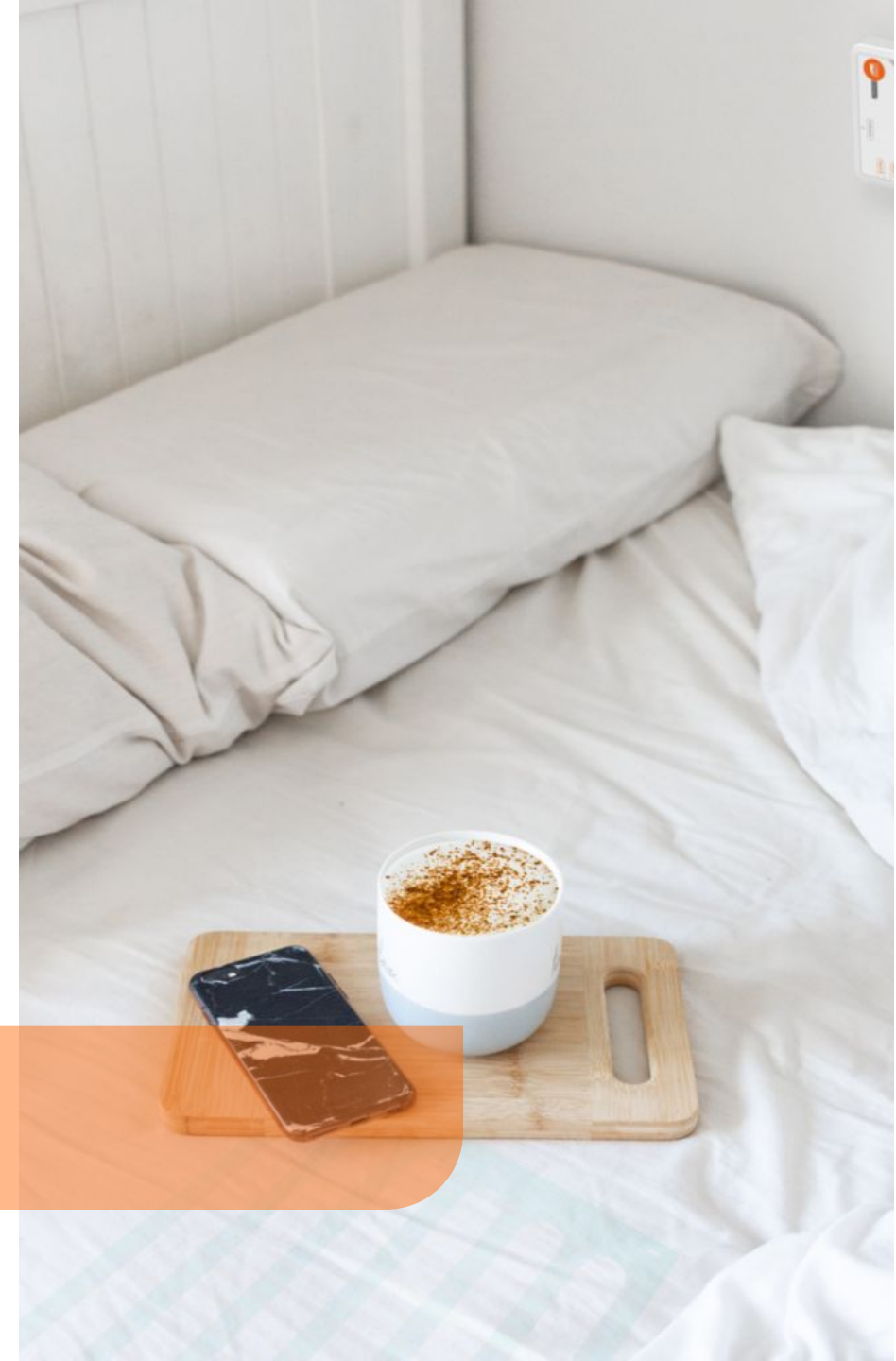


Smart Technology,
Made to Care

Connect

Modes/Options

- No.0. Personal Trigger Mode
- No.1. Bed/Chair Absence Mode
- No.2. Virtual Bed/Chair Mode
- No.3. Property Exit Sensor Mode
- No.4. Virtual Door Mode
- No.5. Pressure Mat Mode
- No.6. Epilepsy Sensor Mode
- No.7. Enuresis Sensor Mode
- No.8. Sound Option - On/Off
- No.9. External Switch Option - Enabled/Disabled (Default - Enabled)

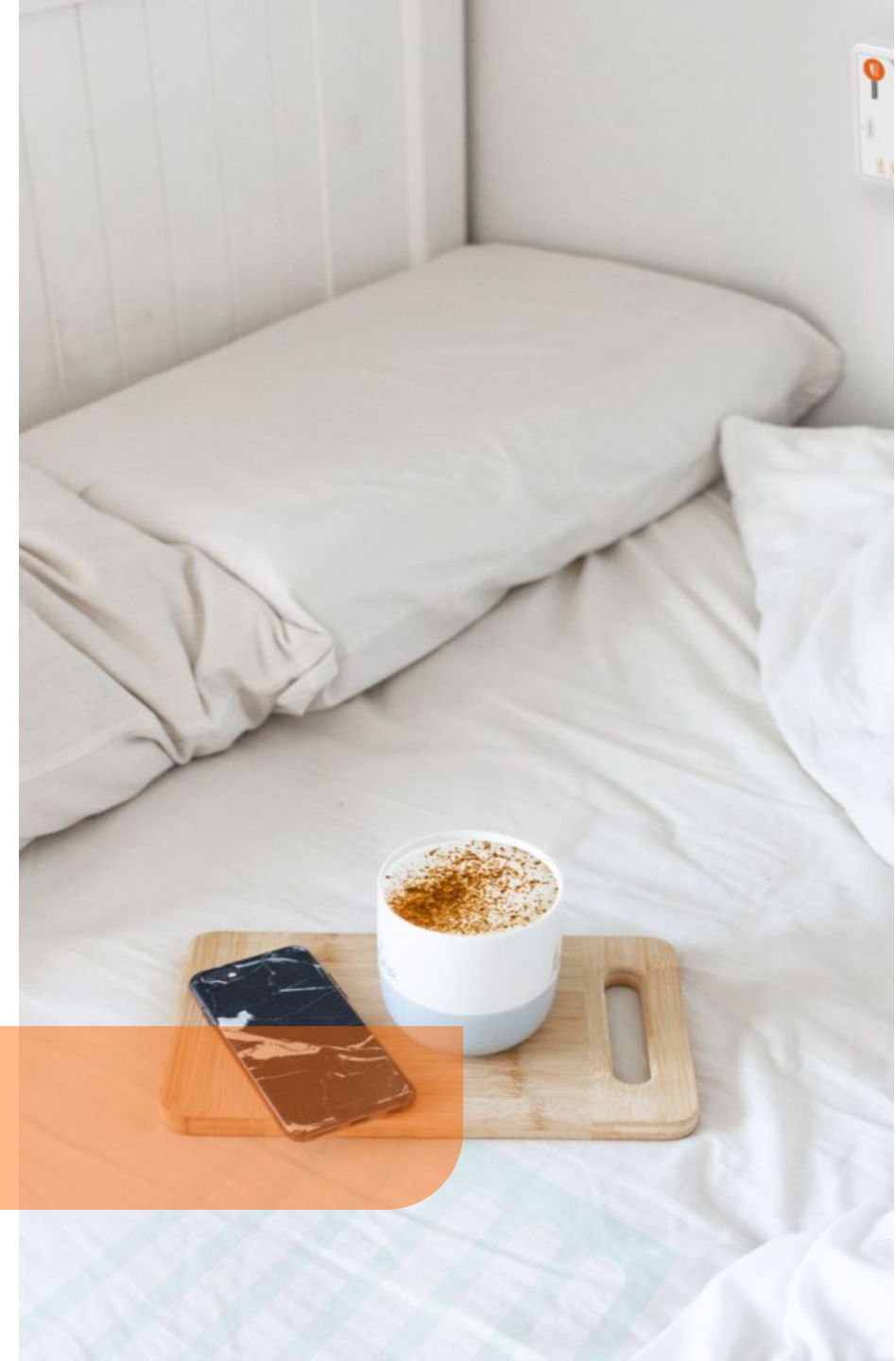


Smart Technology,
Made to Care

Connect

Absence Delay

The absence delay is used for Modes 1-4. Once the modes/options have been chosen, simply turn the dial to the desired number which matches the delay time required. There is no need to press and hold the internal switch afterwards. The time settings are:

- 0 – 0.5 Second
- 1 – 1 Second
- 2 – 3 Seconds
- 3 – 10 Seconds
- 4 – 30 Seconds
- 5 – 1 Minutes
- 6 – 3 Minutes
- 7 – 10 Minutes
- 8 – 15 Minutes
- 9 – 30 Minutes

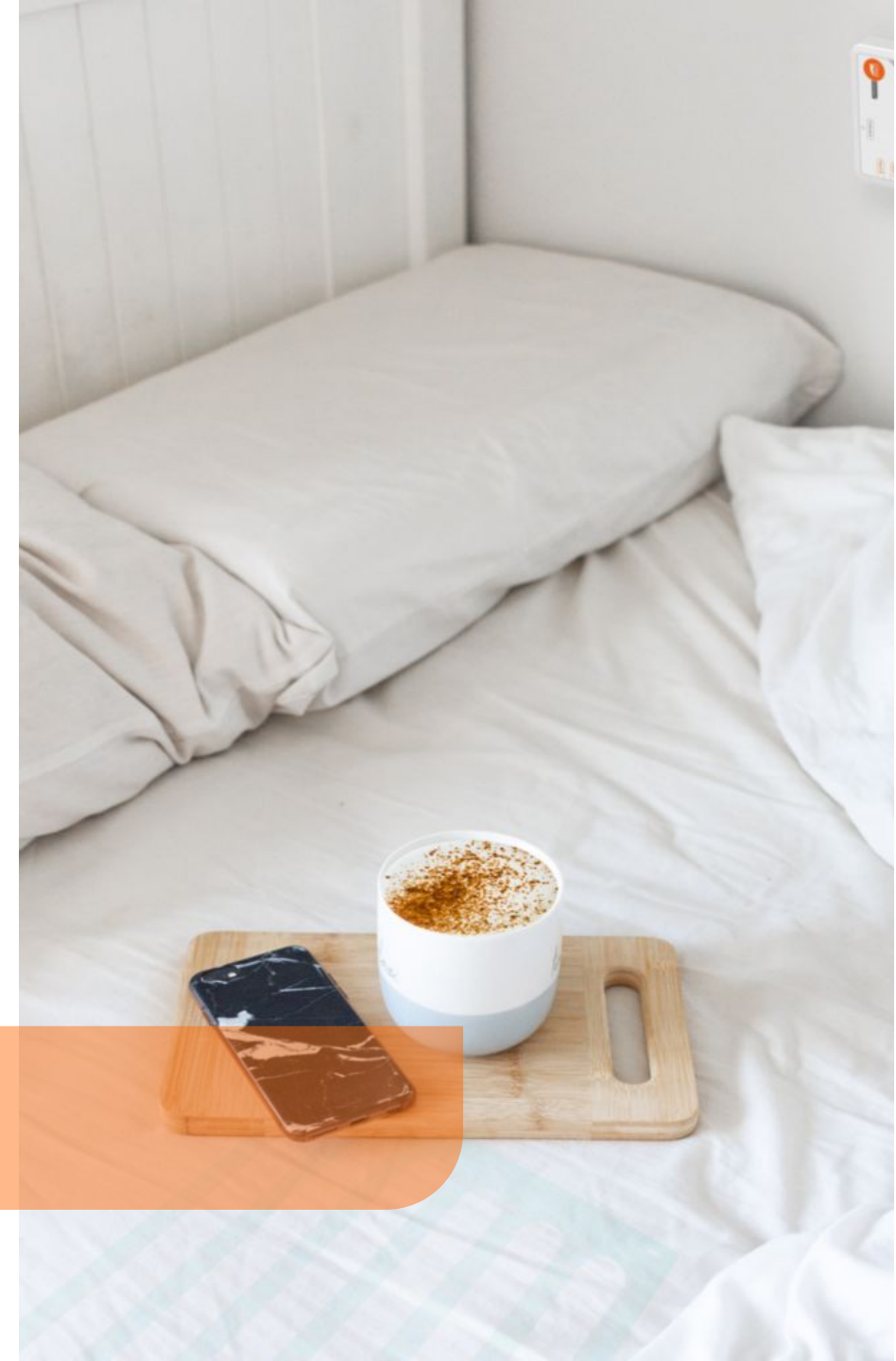


Smart Technology,
Made to Care

[Connect

Registering to a Notifier / Other System

- Program the Connect to the desired Mode/Option (if using modes 1-4, ensure the absence delay is set to a short time).
- Prepare the Notifier/other system by entering registration mode.
- Activate the Connect
- The Notifier/other system should acknowledge the new device.

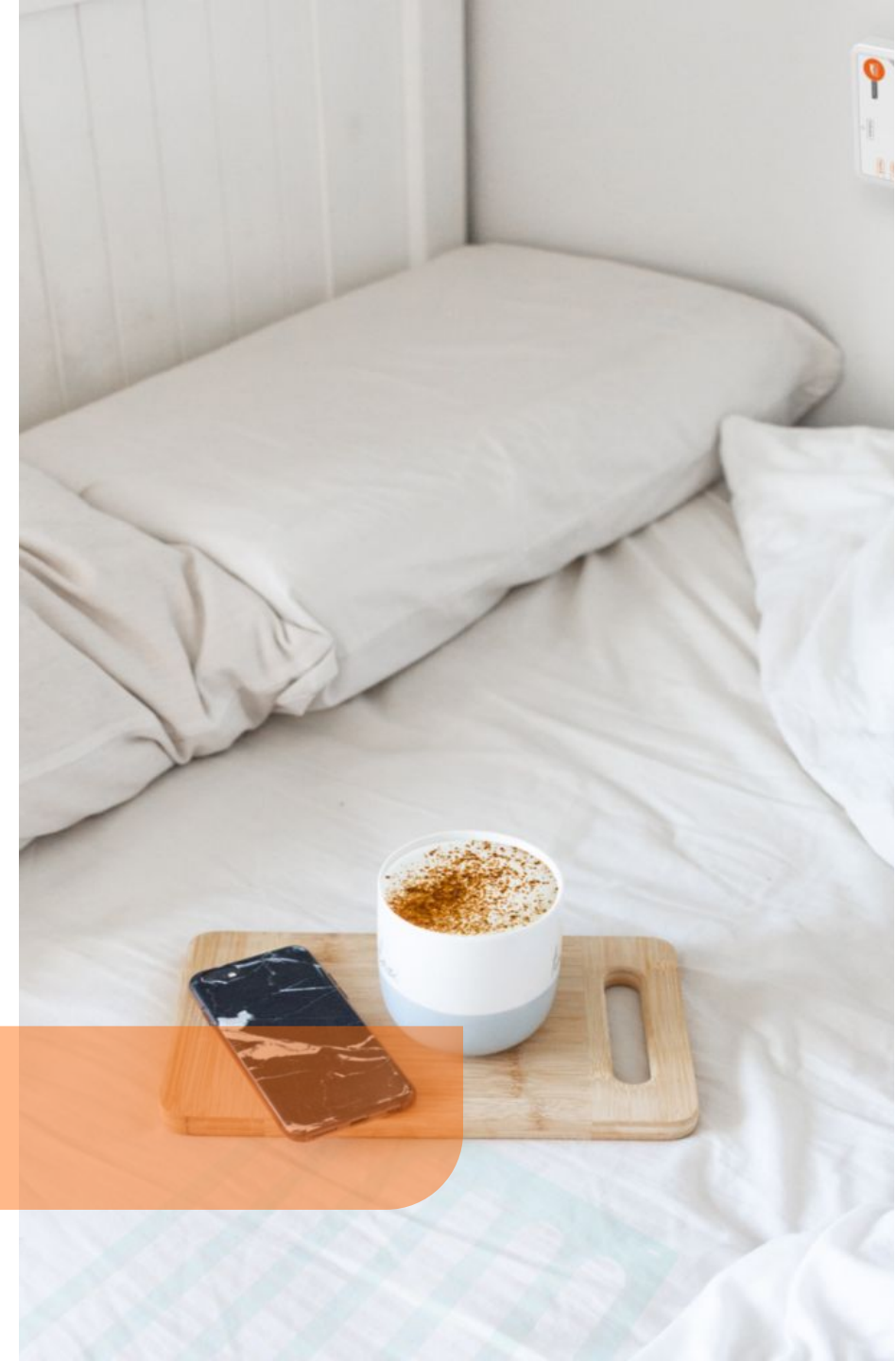

Smart Technology,
Made to Care

Thank you



Smart Technology,
Made to Care

www.sensorium.co.uk

01383 720600 | sales@sensorium.co.uk